



Europäisches Patentamt
European Patent Office
Office européen des brevets



⑪ Publication number:

0 525 577 A2

12

EUROPEAN PATENT APPLICATION

②¹ Application number: 92112402.0

Int. Cl.⁵: **B01L 3/00**, **B01L 9/06**,
G01N 33/53, **G01N 35/02**

② Date of filing: 20.07.92

③ Priority: 26.07.91 US 736155

④3 Date of publication of application:
03.02.93 Bulletin 93/05

④ Designated Contracting States:
DE FR IT

71 Applicant: E.I. DU PONT DE NEMOURS AND
COMPANY
1007 Market Street
Wilmington Delaware 19898(US)

(72) Inventor: **Nurse, Colln A.**
1 Fern Court
Newark, Delaware 19702(US)
 Inventor: **Bernstine, Robert Eric**
111 George Street
Chesapeake City, Maryland 21915(US)


74 Representative: von Kreisler, Alek et al
**Patentanwälte von Kreisler Selting Werner
 Deichmannhaus am Hauptbahnhof
 W-5000 Köln 1(DE)**

⑤4 Carrier device.

57) An analytical carrier device permits the holding of a sample, reaction vessel with reagents for the sample, and a container for receiving a processed

sample. The container facilitates transfer to an analytical device for analysis.

Fig. 3

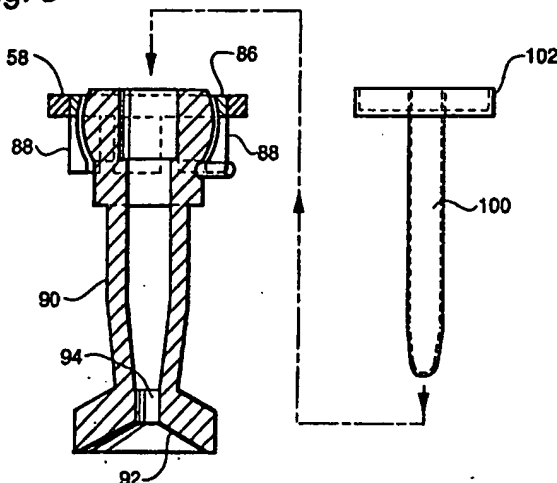

Cross Reference to Related Patent Applications

This invention discloses materials which is disclosed and/or claimed in a patent application entitled Vortex Mixer Drive, Serial No. _____, filed _____ (IP-0801) and also an application entitled Multilinear Automatic Apparatus for Processing Immunoassays, Serial No. _____, filed _____ (IP-0905).

Field of the Invention

This invention relates to a carrier device for holding samples which are to be analyzed in a manner which facilitates their combination with analysis reagents and apparatus for transferring the results for determination by further instrumentation.

Background of the Invention

In the analytical field it is often necessary to process samples of reagents for analysis by combining them with various reagents, support particles and the like. Following such analysis, it is then necessary to transfer the processed sample held in a reaction device back to a device which facilitates processing the results. This is particularly true when the interim processing of a sample involves repeated steps that typically required immunoassay techniques. These include reaction time and wash cycles all involving the use or other particles to facilitate the chemical reactions. In a case of immunoassays, it is necessary to vortex the contents of a reaction vessel to maintain the particles suspended so that the reaction may go to completion. Such analysis may also present difficulties inasmuch as the samples may come from different units and all require mounting on the processing carrier. A further problem arises due to the combination of the samples and reagents during the processing. It is difficult to keep track of the proper sample to make sure it receives the processing needed and does not become mixed up with other samples so that the integrity of the analysis is lost.

Summary of the Invention

Many of these problems are solved by the device of this invention which facilitates the analysis of samples in a reaction vessel. The device comprises a top member having an end portion, a support for the top member, a transparent container having analysis reagents contained therein and having a header containing an instruction code formed on the header, a sample holder being removably mounted on the top member adjacent the header, and a reaction vessel holder mounted by the top member in the end portion, whereby the

sample may be positioned in the reaction vessel for reaction and thereafter transferred to the transparent container for analysis.

In a particularly preferred embodiment, the device of the invention is mounted to permit the sample to be nutated. The reaction vessel itself comprises an inner container having a longitudinal axis and which contains a first reagent and an outer container coaxially positioned about the upper portion of the inner container, the outer container having a second reagent. It is desirable that the transparent container be slidably removeable from the top member to facilitate its use and transfer for further analysis.

The apparatus of this invention maintains three units together, i.e., the sample, the reagent for analysis and the ultimate processed sample in a separate container. This facilitates processing samples and maintains all the units necessary for analysis together during the analysis time interval. Multiple sample cups may be simply clipped on to the carrier for use with different sample input and separate incubation is permitted by the structure of the carrier for material within the reaction vessel.

Brief Description of the Drawings

This invention will be more clearly understood when considered in conjunction with the accompanying drawings in which like reference numerals refer to like components in each of the drawings, in which:

Figure 1 is an exploded view of the carrier device constructed in accordance with this invention;

Figure 2 is a section taken through the stopper of a container along the lines 2-2 of Figure 1; and

Figure 3 is a section taken along the lines 3-3 of Figure 1 particularly depicting the construction of the sample holder.

Detailed Description of the Preferred Embodiment

There may be seen in Figures 1, 2, and 3, exploded sectional views of one of the carriers constructed in accordance with this invention. The carrier is seen to contain a hollow, molded housing 50 defined by a pair of sidewalls 52, a top plate 58, and a base support 60. A drive bar 140 is positioned in the lower portion between the sidewalls and secured to the base support as by glueing. This bar has receptacles 61 to facilitate its receiving driving or positioning pins for positioning the bar 140 and hence the carrier. The housing may be formed of polysulfone or any other suitable engineering plastic which is rigid, strong and chemically inert. Attached to the front sidewall (in the

drawing) is a partition 54 which cooperates with the top 58 to accommodate the top frame of an analytical pack 62 of an analytical pack 64 which may be the same and preferably is the same as the aca® pack used in the aca® Automatic Clinical Analyzer sold by E. I. du Pont de Nemours and Company, Wilmington, Delaware, U.S.A. The aca® pack has identifying indicia 66 on the top which may be read by appropriate sensors to indicate the particular test being run and includes a hollow septum 68 with an orifice 70 which may be used to introduce materials into the a plastic pack 72. Since the aca® pack is well known it will not be described further.

In any event, the partition 54 and top 58 cooperate to define an orifice 56 adapted to accommodate the top member of the aca® pack 62 so it may be inserted into the carrier with the lower side pack 72, which is formed of plastic material. The side pack is to slide in between the two walls 52. The top of the carrier 50 also includes an elongated cup-like member 76 which is adapted to receive a removable sample reservoir 78 containing a reservoir 80. The sample reservoir 78 is held in the position within the opening 76 by appropriate molded grips 82. A fitting feature 84 may be provided for the sample holder 78 to control access to the opening.

To complete the carrier 50, the end of the top member 58 may have an orifice 86 with downwardly extending flanges 88 adapted to hold a reaction vessel holder 90. The flanges 88 are concave on the inside to define a socket which cooperates with the bulbous top on a reaction vessel holder 90 in a ball and socket joint manner. The lower portion of the reaction vessel holder 90 may be shaped as to have an inverted cavity or receptacle 92 at the upper end of which is a bore 94 adapted to receive a pin from a nutating drive member.

In an alternative embodiment of this invention, the reaction vessel holder 90 may be the reaction vessel itself although the use of the holder is preferred for its long term stability and reliability. If the reaction vessel 90, as a tube holder, is adapted to receive a reaction vessel 100, the vessel has at the upper portion thereof a concentric chamber 102 for holding reaction reagents that typically may be used, for example, in an immunoassay process.

The reaction vessel holder 90 is nutated by an automatic apparatus 104. Although any suitable drive apparatus may be used that provides two directions of linear motion and one direction of rotary motion, that described in copending application Serial No. _____, filed _____, (IP-0905) is one that has been used successfully. This drive apparatus provides a bidirectional motion as depicted by the line 106 (Fig. 1) as well as rotational motion as depicted by the line 108 to the reaction vessel holder 90. The drive apparatus is

powered by a single bidirectional drive motor 110 which provides rotational motion to the drive apparatus 104. The automatic apparatus engages the reaction vessel holder 90 by elevating a mixing cylinder on which a pin is positioned contiguous the periphery at a point off the elongated axis of the mixing cylinder. In other words the pin engages the bottom end of the mixing vessel 90 in a position which is eccentric to the axis which mounts the mixing cylinder. The apparatus then spins the cylinder moving the engaged end of the vessel into an orbit. If the vessel is managed so that it is free in two rotational directions of freedom, then the contents of the reaction vessel holder 90 will swirl or nutate thus mixing them. Reversal of the drive which spins the mixing cylinder 110 stops the orbiting of the vessel and lowers the cylinder thus disengaging the cylinder from the reaction vessel holder 90.

The carrier device just described has many advantages and it permits in effect the combination of three different functions into a single device. Firstly, sample cups from various analysis devices holding the sample may be snapped into place in the upper portion of the device. Secondly, a reaction tube containing reagents for the analysis of samples may be inserted into a carrier tube which may be inserted into and processed separately from the carrier unit. Thirdly, the completed reacted sample with reagents may be introduced into a separate unit which is slidably removeable from the carrier for subsequent processing and analysis is desired in other machines.

Claims

1. An analytical carrier device for facilitating the analysis of samples in a reaction vessel comprising:
 - a top member having an end portion,
 - a support for the top member,
 - a transparent container having analysis reagents contained therein and having a header containing an instruction code formed on the header,
 - a sample holder being removably mounted on the top member adjacent the header, and
 - a reaction vessel holder mounted by the top member in the end portion, whereby the sample may be positioned in the reaction vessel for reaction and thereafter transferred to the transparent container for analysis.
2. An analytical carrier device as set forth in Claim 1 wherein the sample holder is mounted rotatably to permit its lower end to be nutated.
3. An analytical carrier device as set forth in

Claim 1 wherein the reaction vessel holder includes a reaction vessel positioned in the holder.

4. An analytical carrier device as set forth in Claim 1 wherein the reaction vessel comprises an inner container having a longitudinal axis and which contains a first reagent and an outer container coaxially positioned about the upper portion of the inner container, the outer container having a second reagent.

5

10
5. An analytical carrier device as set forth in Claim 1 wherein the reaction vessel holder includes a reaction vessel partially rotatable in the holder.

15
6. An analytical carrier device as set forth in Claim 1 wherein the reaction vessel comprises an inner container having a longitudinal axis and contains a first reagent and an outer container coaxially positioned about the upper portion of the inner container, the outer container having a second reagent.

20

25
7. An analytical carrier device as set forth in Claim 1 wherein the transparent container is slidably removable from the top member.

30
8. An analytical carrier device as set forth in Claim 1 wherein the transparent container is slidably removable from the top member.

35

40

45

50

55

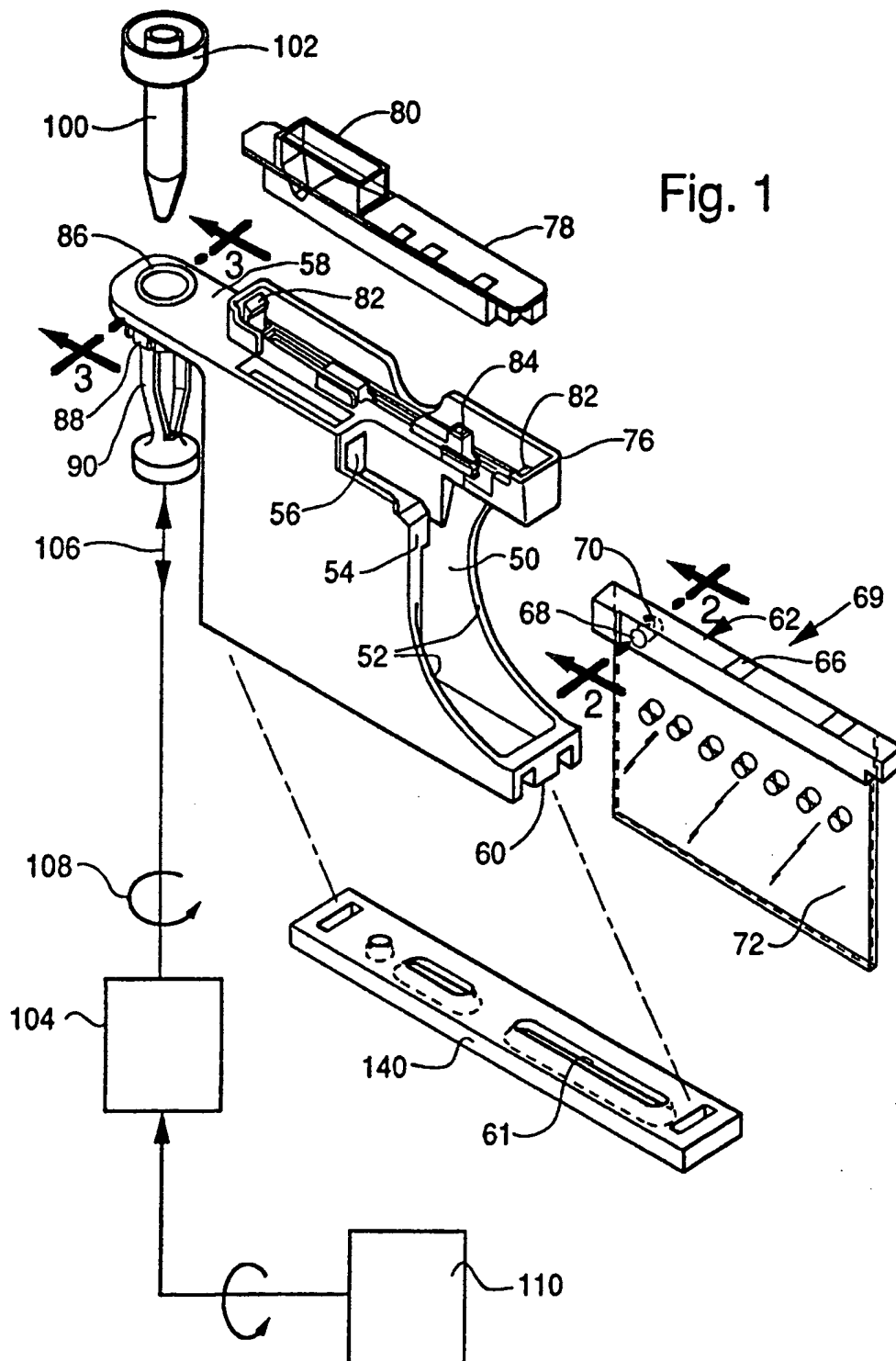


Fig. 2

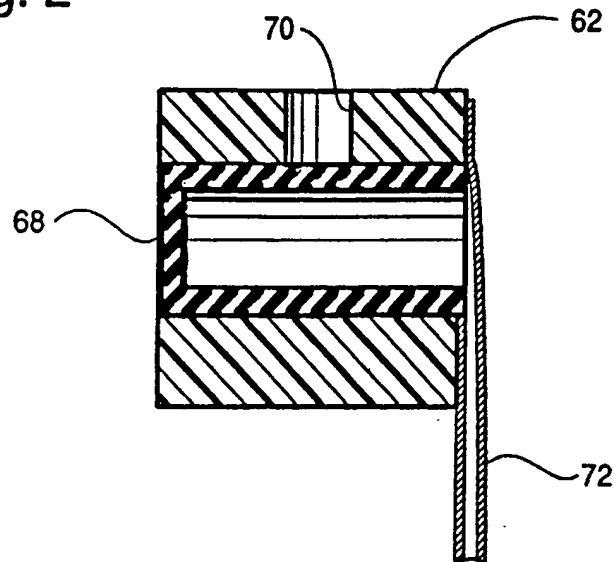


Fig. 3

